


展开 520*150mm 折后130*150mm
公差 ±1mm

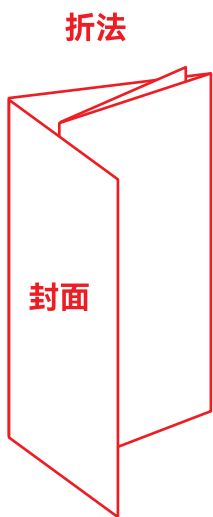

封底

封面

08

11-13mm 13-15mm 15-17mm 17-20mm 20-23mm 23-26mm

15mm 17mm 19mm 21mm 24mm 27mm

09

Placement in Bra

① Place the breast shield inside your bra, ensuring the nipple is centered.

② Secure the pump in place for hands-free operation.

Ensure that the breast pump fits snugly against the breast without any gaps.

EZMOM

www.ezmom.com

EZMOM

Wearable Breast Pump

EZ-NATURE

User Guide

Operating Guide

Nipple Measuring Guides

If you have any questions, please contact us by email.
support@ezmom.com

01

A 24mm Double-sealed Flange B Membrane C Valve D Milk collector (180ml) E Pump F Dust cover

02

a: On/Off button / Mode button
b: Increase vacuum
c: Reduce vacuum

- “ ” = Massage mode
- “ ” = Express mode
- “ ” + “ ” = Mixed mode

03

04

Charging Indicator:

purple breathing light	1%-21%
orange breathing light	22%-50%
green breathing light	50%-99%
a steady green light when fully charged.	

Remaining Battery Indicator:

The battery level is displayed on the battery icon.

purple light flashes	1%-3.5%
purple light	3.5%-22%
orange light	22%-50%
green light	50%-100%

05

06

Replace washable parts 1-3 months

07

Keep parts away from the heated container's edges.

08

Nipple Measuring (Don't skip this step)

A properly fitted flange maximizes milk output, prevents nipple damage, and ensures effective milk pumping.